



# XenoSure® Biologic Vascular Patch

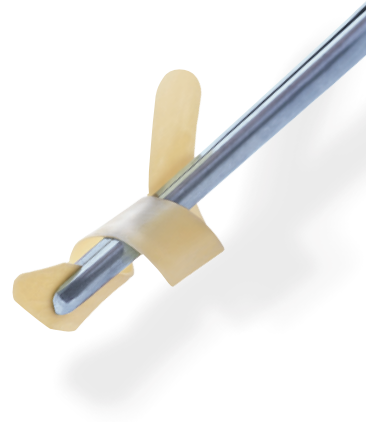

# XenoSure® Biologic Vascular Patch

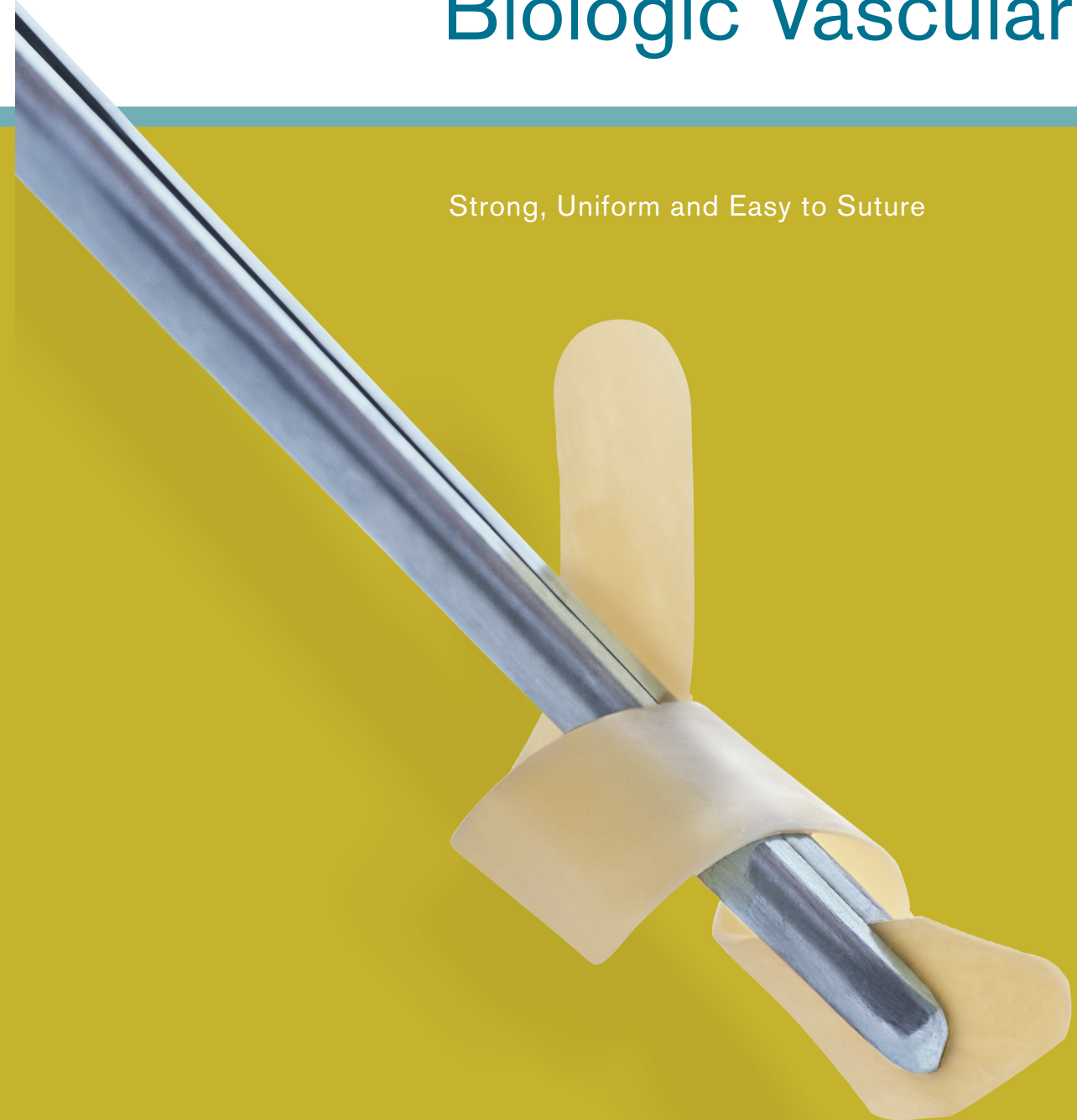

Strong, Uniform and Easy to Suture

## ORDERING INFORMATION

	Width	Length	Model #
<b>XenoSure Biologic Vascular Patch (Bovine Pericardium)</b>			
Thickness range 0.350–0.750mm			
XenoSure Biologic Vascular Patch	1 cm	6 cm	e1P6
	1 cm	14 cm	e1P14
	0.8 cm	8 cm	e0.8P8
	2 cm	9 cm	e2P9
<b>XenoSure Biologic Surgical Patch (Bovine Pericardium)</b>			
Thickness range 0.350–0.750mm			
XenoSure Biologic Surgical Patch	4 cm	4 cm	e4P4
	4 cm	6 cm	e4P6
	6 cm	8 cm	e6P8
	8 cm	14 cm	e8P14
	10 cm	16 cm	e10P16

These specifications are not intended as a warranty. In the interest of product improvement, these specifications may be changed from time to time without notice. Please consult your sales representative for details.

### LeMaitre Vascular, Inc.

63 Second Avenue  
Burlington, MA 01803  
USA

T 781 221 2266  
F 781 221 2223

### LeMaitre Vascular ULC

5915 Airport Road, Suite 608  
Mississauga, ON L4V 1T1  
Canada

T 905 673 2266  
F 905 673 2223

### LeMaitre Vascular GmbH

Otto-Volger-Str. 5a/b  
65843 Sulzbach/Ts.  
Germany

T +49 (0)6196 659230  
F +49 (0)6196 527072

### LeMaitre Vascular GK

1F Kubodera Twin Tower Bldg.  
2-9-4 Kudan-minami, Chiyoda-ku  
Tokyo 102-0074  
Japan

T +81 (0)3 5215 5681  
F +81 (0)3 5215 5682

### Customer Service

US 800 628 9470  
Canada 855 673 2266

[www.lemaitre.com](http://www.lemaitre.com)



LeMaitre, XenoSure and Your Peripheral Vision are registered trademarks of LeMaitre Vascular, Inc. ©2015 LeMaitre Vascular, Inc. All rights reserved. Printed in the USA. M0236 Rev. J 04/15



Your Peripheral Vision®

# XenoSure® Biologic Vascular Patch

## STRONG, UNIFORM AND EASY TO SUTURE

The XenoSure Biologic Patch is a high-quality bovine pericardium patch used for precision endarterectomy and vascular reconstruction. Using the same tissue technology developed for heart valves, the XenoSure patch is strong, uniform and easy to handle and suture.



### BENEFITS:

- Biocompatible
- Excellent tensile and suture retention strength
- Does not require special sutures
- Uniform collagen thickness\* for easier suturing

### APPLICATIONS:




- Carotid endarterectomy
- Femoral, iliac, renal and tibial patching
- Profundaplasty
- Arteriovenous access revisions

### HANDLING:

- Similar to autologous tissue
- Excellent handling characteristics
- Easily trimmed to desired shape or size

## Available Sizes

### XENOSURE BIOLOGIC VASCULAR PATCH

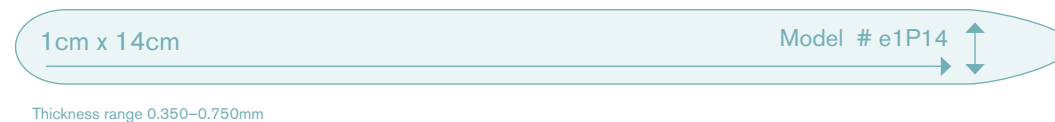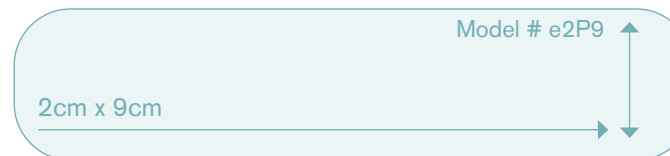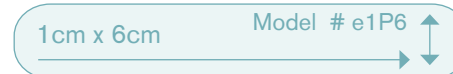


### Rinse Procedures 4cm x 6cm & Smaller

1. Fill basin with 500ml sterile saline
2. Remove with sterile atraumatic forceps
3. Agitate patch in saline for 2 minutes
4. Leave in saline until required by surgeon

### Rinse Procedure 6cm x 8cm & Larger

1. Fill 1 basin with 1000ml sterile saline
2. Remove with sterile atraumatic forceps
3. Agitate patch in saline for 3 minutes
4. Leave in saline until required by surgeon

### XENOSURE BIOLOGIC SURGICAL PATCH

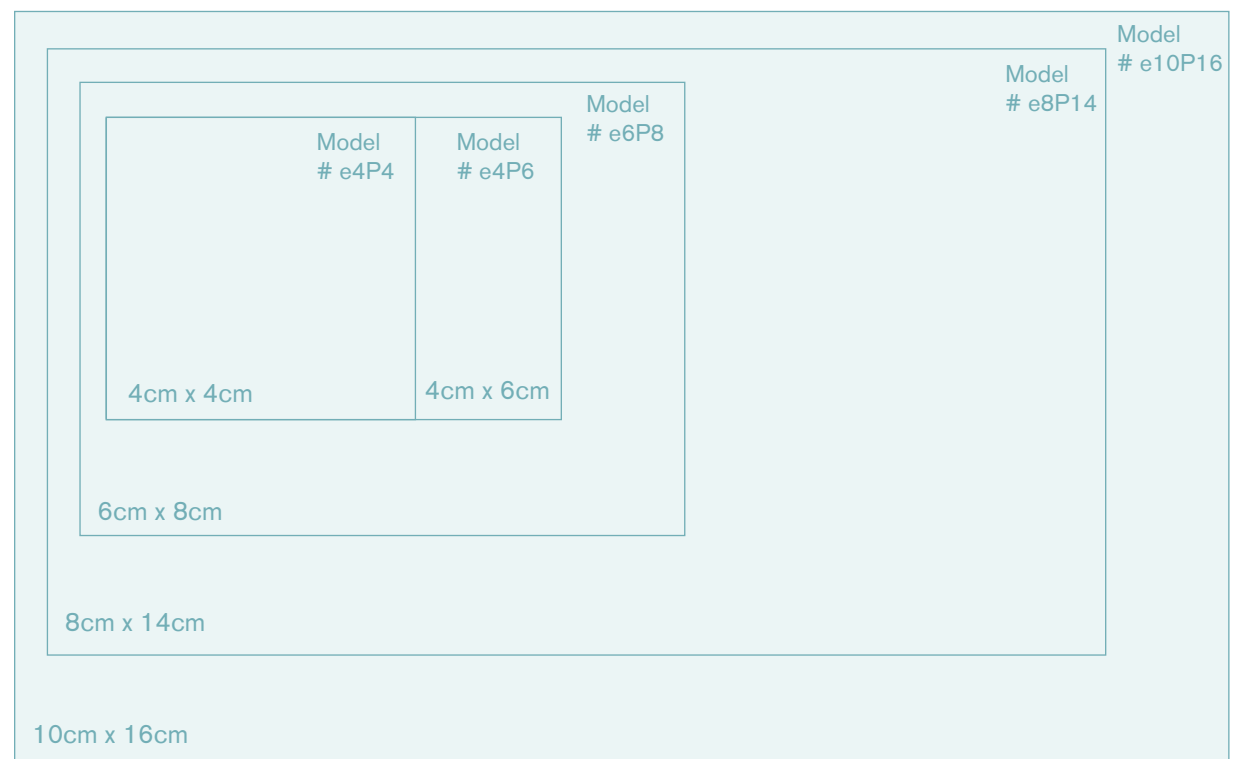

Thickness range 0.350-0.750mm

\* Data on file at LeMaitre Vascular, Inc.